


# Sample Preparation



## Grids

**Thickness:** 25  $\mu\text{m}$   
**Diameter:** 3 mm  
Film on “dull” side  
“Shiny” side has the metal only

- Ultrathin Carbon Film on a Lacey Carbon Support Film
- [https://www.tedpella.com/Support\\_Films\\_html/Support\\_Films\\_and\\_Substrates\\_Overview.aspx](https://www.tedpella.com/Support_Films_html/Support_Films_and_Substrates_Overview.aspx)
- 01824 Ultrathin Carbon Film on Lacey Carbon Support Film, 400 mesh, Copper
- 160 PELCO® TEM Grid Storage Box

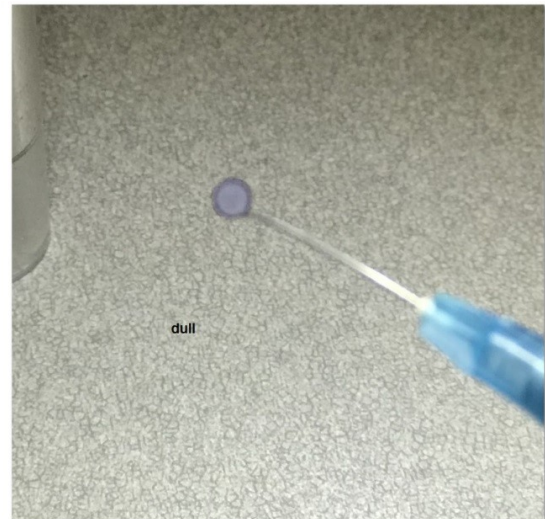


## Basic Shake and Bake

Put a small amount of powder in Eppendorf tub (enough to coat the bottom). Put one grid into the tube. Close the tube and shake for 1 minute. Open the tube and remove the grid. The grid should have enough material on it and is ready to be imaged in the TEM.

## Instructions for Nanoparticle Dropcast TEM sample preparation

### Things you will need:

- Ethanol (or water, IPA, or dimethyl carbonate)

- A small metal scooper
  - A sonicator
  - Several small vials that are centrifuge compatible
  - A centrifuge
  - Nanoparticles
  - UltraThin 3 nm Carbon Grids
  - Glass slide
1. Check the MSDS for mixture with ethanol (make sure the nanoparticles are NOT pyrophoric)
  2. Use the small metal scooper to extract about ½-1 cms worth of nanoparticles
  3. Place the nanoparticles in a small vial
  4. Fill the small vial with ethanol up to the 1 cm below the rim
  5. Sonicate the filled small vial for 10 minutes on high power (make sure to find an anti-node)
  6. Centrifuge for 20 minute (make sure to balance the centrifuge)
  7. Use a pipette to collect several uL of fluid just above the concentrated part at the bottom of the vial. See the arrow below . The pipette should just touch above dark part
- 
8. Deposit 1 uL of fluid onto the dull (not shiny) side of the grid
  9. Let the grid dry on a glass slide
  10. Vacuum bake at 50-400 degree C for several hours if contamination is found of the surface
  11. Plasma clean as needed

From:  
<https://bpm-wiki.cnsi.ucsb.edu/> - NSF BioPACIFIC MIP Wiki

Permanent link:  
<https://bpm-wiki.cnsi.ucsb.edu/doku.php?id=microed-sample&rev=1727905305>

Last update: **2024/10/02 21:41**

