2025/11/06 03:52 1/9 EinScan-SP Training SOP

EinScan-SP Training SOP

Tool Name and Location

The SHINING3D EinScan-SP 3D Scanner is located in Fume Hood #1 in Elings 2448.

Safety

Machine Limitations:

 The turntable has a maximum weight limit of 5kg (~11lbs). Exceeding this weight could damage the machine and its internal components.

Environmental:'

• Keep the scanner in a well ventilated area. Keep flammables, corrosive agents, and liquids away from the scanner. Failing to do so will damage the machine and its internal components.

Physical:

• The projector emits a bright light. Avoid looking directly into the light source to prevent eye injury.

Training

training requirements: Users must attend the in-person training session and should consult the Safe Operating Procedure and EinScan-SP Manual (both found on the Innovation Workshop GauchoSpace page). **Training Outline**

- 1. Overview of the scanner and the system setup
 - 1. Structured light projector and dual camera setup
 - 2. Electronically controlled turntable
 - 3. Laptop computer
- 2. Safety
 - 1. Machine Limit
 - 1. Do not exceed 5 kg (11 lbs)
 - 2. Hazards
 - 1. Keep well ventilated, away from corrosives, liquids, and flammables.
- 3. Walkthrough of starting up EinScan-S software
 - 1. File formats, different scan types
 - 2. EinScan-S menu
- 4. Scanner Limitations (0.03 mm resolution)
 - 1. Bigger than 30x30x30 mm, smaller than 250x250x250 mm, less than 5 kg
 - 2. Do not scan reflective, transparent, or fuzzy objects
 - 3. Paint or coat dark objects with a light color
- 5. Step by step walkthrough of a scan
 - 1. Scanning with the turntable
 - 2. Scanning without the turntable

- 3. Adjusting the brightness and stitching scans
- 6. Software walkthrough
 - 1. Turntable toolbar
 - 2. Scan editing toolbar
 - 3. Main toolbar (alignment, scans, exporting, brightness)
- 7. Editing, Meshing, and exporting completed scans


Safe Operation of this Tool

Before Using the Scanner: Make sure that the object:

- Is larger than 30x30x30mm (1.2×1.2×1.2in)
- Is smaller than 250x250x250mm (38x38x38in)
- Weighs less than 5kg (11lbs)
- Isn't shiny, transparent, or furry
- · Has all dark features lightened with a powder or coating

Sequence of Operation:

- 1. Make sure the scanner is powered on and open the EinScan-S software (on desktop)
- 2. Select EinScan-SP
- 3. Under "Choose Working Mode", select Fixed Scan
- 4. Select <u>New Scan</u> or <u>Open Scan</u>. Name the file and save to your group's folder to define the location for saved scan images (saved as .fix prj file)
- 5. Select Textured Scan or Non-Textured Scan
 - 1. Textured scan is used to scan colors and fine features. The user must run a white balance test when using this feature.
- 6. Adjust <u>Brightness slider</u> until no red areas are visible on the desired part

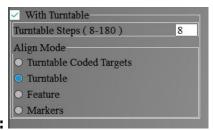

Scanning With the Turntable:

- 1. Check the <u>Use Turntable</u> box and specify the number of increments (8-180 steps). The <u>Align Mode</u> menu will appear.
- 2. Under Align Mode, select the desired alignment method
 - 1. Turntable
 - 2. Features
 - 3. Markers (External markers must be placed on the surface of scanned object)

Scanning Without the Turntable:

1. Place the object in the desired orientation within frame of the two cameras

Starting and Editing Scans:

2025/11/06 03:52 3/9 EinScan-SP Training SOP

- 1. Select <u>Start Scan</u> . Press <u>Pause</u> or <u>Delete Current Scan Data</u> if necessry
- 2. Use Shift + Left Click and the toolbar to remove any unwanted scan data. Press the <u>Green Check</u> to confirm.
- 3. Reorient the part and continue to scan until object is fully modeled (Same controls to modify the stitched model)
- 4. Correct misaligned stitches by selecting <u>Align</u> and using Shift+Left Click to choose three common points to manually correct the software.

Exporting and Saving the Scan:

- 1. When the object is fully scanned, select Mesh
- 2. Select <u>Watertight</u> to create a closed object file or an <u>Unwatertight</u> model to directly edit the object file
- 3. After the object has been meshed, specify the desired <u>Simplification Ratio</u>. Make sure that <u>Smooth</u> and <u>Sharpen</u> are both selected. Specify <u>Hole Fill Size</u> if creating an Unwatertight model
- 4. Select <u>Save</u> and specify a <u>Scaling Ratio</u>. Save the file as a (.obj .stl .ply .3mf or .asc if you wish to edit the saved mesh)

Editing the Exported File: We do not currently have a universal software solution to edit and modify the generated part files. Please contact the lab staff if you have a solution to this problem. We will continue to research solutions. Meshmixer by Autodesk is so far the most promising mesh file editor, but is not available on lab computers for now.

Detailed Procedure-Step by Step

Supplies:

No supplies are necessary to use the EinScan-SP 3D scanner. A blackout enclosure is supplied to prevent light from affecting the scan quality.

Setup:

The default Innovation Workshop setup is to have both the scanning unit and the turntable attached to the scanner stand. Please contact the staff if you wish to try scanning from the tripod instead. The scanning unit and the turntable are both attached to the scanner stand.

The Scanning Unit:

Contains two cameras and a projector that is capable of projecting both structured and RGB light. The offset cameras provide two perspectives of the same point, allowing the software to calculate depth by comparing the two images. The cameras also capture the distortion of the structured light vertical lines projected over the object surface. This similarly measures the depth and features of the surface. The cameras also superimpose the images taken with Red, Green, and Blue light to create a color accurate 3D model of the object. This color image is combined with the structured light depth scans to create a colored 3D model. The scanning unit can also be removed from the stand and attached to the provided tripod to scan larger objects. The turntable is connected to the scanning unit via USB cable, and the scanning unit is connected to the laptop computer running EinScan-S software via USB cable. The power cable plugged into the scanning unit provides power for both the scanner and the turntable.


The Turntable:

Incrementally rotates 360 degrees per scan to capture a complete image of the object. The degrees rotated during each increment are adjustable in the EinScan-S software and range from 2 to 45 degrees. The turntable is controlled and powered via USB connection to the scanner. The surface of the turntable contains coded markers, which are used to reference the position of the object when the scanner is detached from the stand and placed on the provided tripod.

Operation: Prepare the object for scanning. Use the following checklist


- The part is larger than 30x30x30mm (1.2×1.2×1.2in) and smaller than 250x250x250mm (38x38x38in)
- The part weighs less than 5 kg (11 lbs)
- Make sure there are no transparent surfaces, reflective surfaces, or finely textured surfaces (like hair) on the object to be scanned. These surfaces affect the structured light and decrease the quality and resolution of the scan.
- Powder or coat any dark surfaces. This well help highlight the dark surfaces and increase the quality of the scan.
- Check that the scanner is powered on. The power button is located on the back of the scanner.

 A glowing blue border outlines the power button when the scanner is turned on
- From the desktop, launch the EinScan-S software . All scanning functions are controlled from this software, and will not initialize properly if the scanner is not powered on.

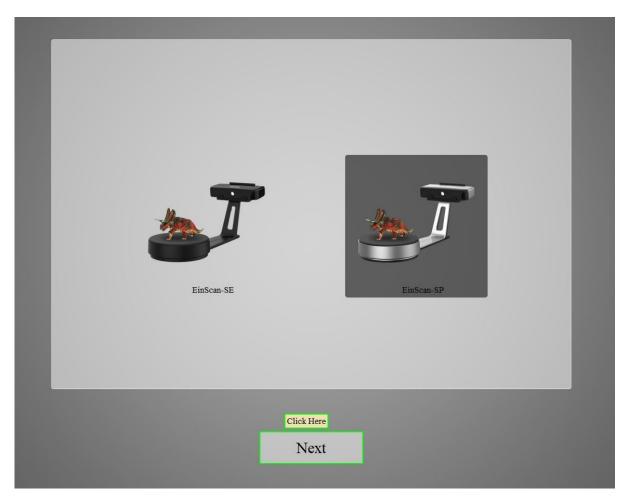


• Once the software has initialized, the scanner menu will pop up. Select EinScan-SP, and then click Next.

2025/11/06 03:52 5/9 EinScan-SP Training SOP

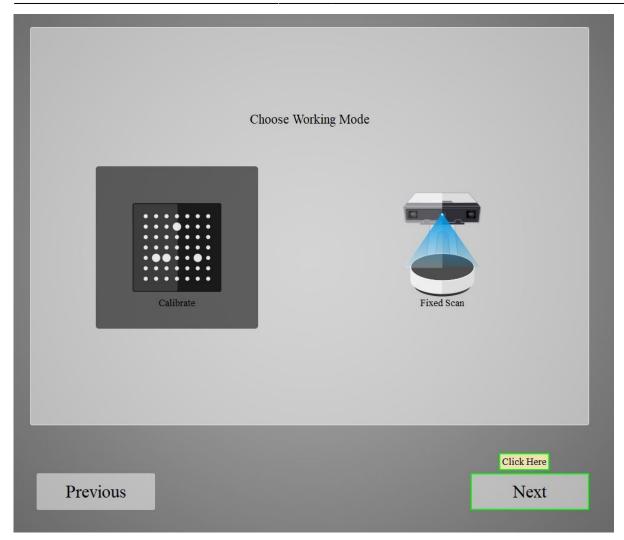

• The working mode menu will pop up. Select Fixed Scan and then click Next

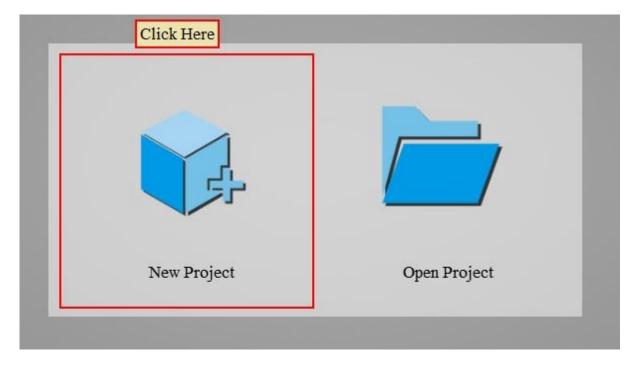


- The projector will start emitting white light to shine on the part.
- Select New Project or Open Project and click to initiate.

To Start A New Scan:

• Name the project and designate a location. This will create a folder to store the scan images and the stitched project (saved as a .fix_prj) file. Individual images can be selected and deleted from this folder


To Open a Scan:

• Locate the project folder and open the .fix_prj file to open the scan project. The file will rebuild the scan using the images present in the folder.

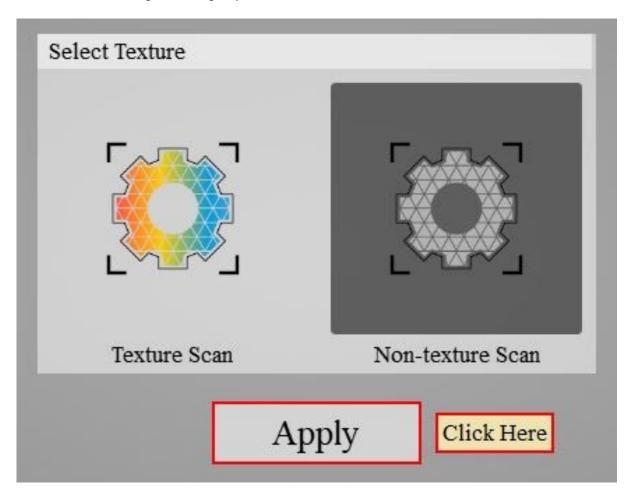



- The texture menu will appear, and give the option of either <u>Texture Scan</u> or <u>Non-texture Scan</u>.
 - A Texture Scan projects colored light in addition to structured light to create a color accurate scanned model of the object. Due to the additional images taken and data processing that must occur, the textured scan process is significantly longer. This mode is not recommended if you simply want to generate a 3D object.
 - A Non-texture scan only uses structured light to capture the surface of the object to be scanned. No information about the color of the surface is captured, although colors on surface may show up separately due to the differing contrast (ie light colored text on a darker background will show up blotchy, even if the surface is flat). The non-textured scan process takes significantly less time than the textured scan process, and is recommended for scans where surface color is not important.

2025/11/06 03:52 7/9 EinScan-SP Training SOP

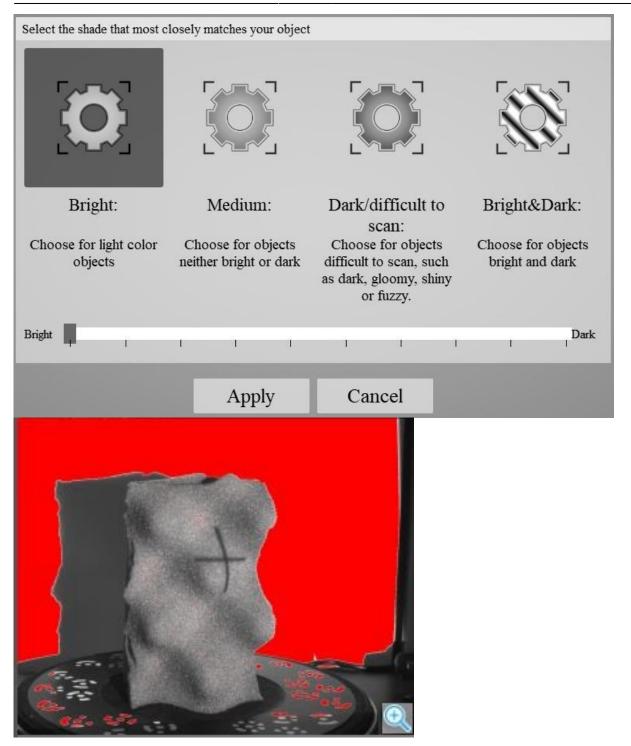

Once a texture mode is selected, the brightness slider will appear. Use the camera screen in the
upper right as a reference for adjusting the brightness level. Move the slider so that the
crosshair is clearly illuminated on the object in the view window, but not so brightly illuminated
that red areas show up on the object to be scanned. An example image of the right exposure
over the part is shown.

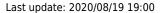

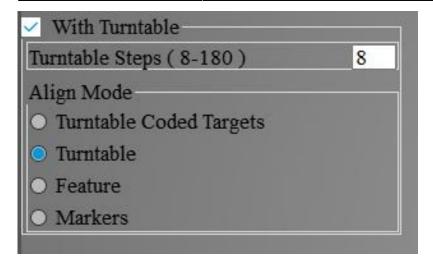

- Scanning with turntable
 - Make sure that the option With Turntable is checked. The following instructions are written under the assumption that the system is mounted on the scanning stand.
 - Make sure that the option With Turntable is checked. The following instructions are written under the assumption that the system is mounted on the scanning stand.
- Tripod Scanning Methods
 - Turntable Coded Targets, Feature, and Markers all require detaching the scanning unit from the scanner stand and mounting it on the provided Shining3D tripod. Details on these methods and using the tripod scanning setup are found HERE.
- Scanning without Turntable
 - Make sure that the With Turntable box is unchecked. The Align Mode toolbar will
 minimize. This mode is used for conducting single image scans, instead of a complete
 rotation of the turntable. This mode is helpful for capturing specific features, or if only one
 surface of an object is desired to be scanned.
- Once the desired alignment mode is set and the brightness has been properly adjusted. It is time to start the scan (Note: The use of the turntable, the number of turntable steps, and the brightness can all be changed from scan to scan. This is a useful way to capture different features and create a more complete scan of the object.)
- The scan is started by clicking Start Scan . This icon is also used to restart a paused scan.

The scan can be paused at any point by clicking Pause Scan . This will halt the scan, allowing you to continue or delete the current scan data. Clicking Delete Scan Data will abort the current scan and remove any captured data. If viewing the stitched object with multiple scans, selecting this icon will delete the most recent scan.

Accident Response:

From:

https://bpm-wiki.cnsi.ucsb.edu/ - NSF BioPACIFIC MIP Wiki

Permanent link:

https://bpm-wiki.cnsi.ucsb.edu/doku.php?id=einscantrainginsop&rev=1597863608

Last update: 2020/08/19 19:00



2025/11/06 03:52 9/9 EinScan-SP Training SOP