


LumiDox II LED Array

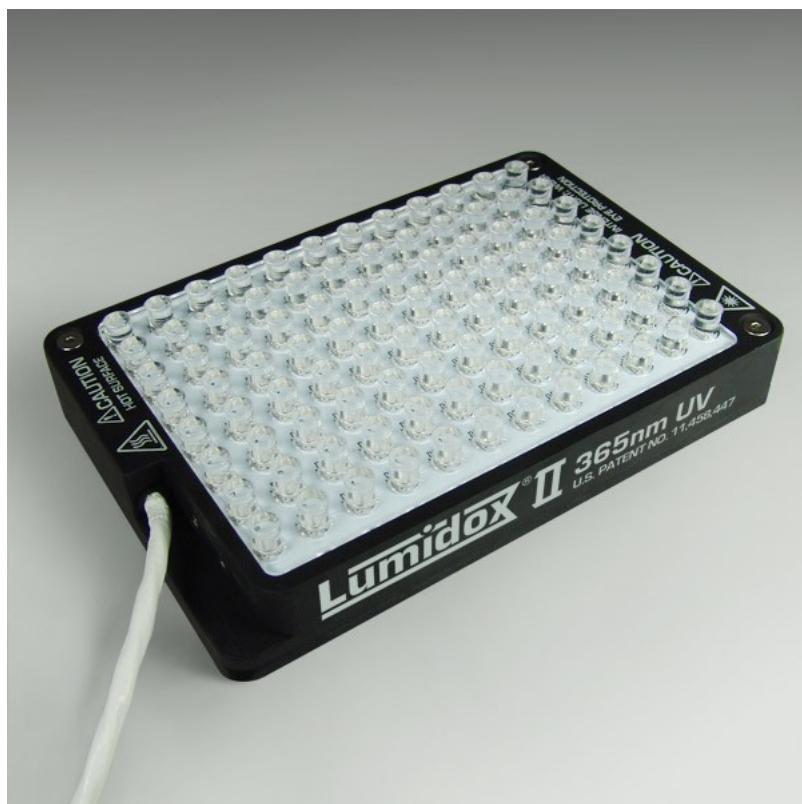

Tool Type:	96-well LED Array
Location:	Elings Hall 2411
Manufacturer:	Analytical Sales
Description:	Chemspeed and Zeiss Microscope Compatible
Principal Scientists	E-mail
Juan Manuel Urueña	jmuruena@ucsb.edu
Morgan Bates	morganbates@ucsb.edu

About

The LumiDox Gen II LED Plate is an innovative tool designed for precise light delivery to 96-well plates. This system utilizes high-intensity LEDs arranged to provide uniform, controlled light exposure across all wells, making it ideal for light-activated experiments. The plate allows for customizable light intensities and wavelengths, supporting a variety of experimental conditions. The LumiDox Gen II's user-friendly interface and robust design enable consistent and repeatable illumination, advancing research in areas such as photochemistry, photobiology, and cell signaling.

Detailed Specifications

- * LED technology:** High-intensity LEDs
- * Compatible formats:** 96-well plates and 96-well shell vial plates (1 mL scale)
- * Wavelengths available:** 365, 420, 445, and 530 nm

* **Power control:** Adjustable power settings

* **Uniformity:** High across all wells

Safety Concerns

The LumiDox Gen II LED Plate produces high-intensity light, which requires the following safety precautions:

- Always inspect the LEDs before use to ensure no damage or malfunction.
- Do not stare directly at the light sources during operation.
- Use protective eyewear rated for specific wavelengths, especially when working with high-intensity light in the UV range.
- Ensure that all users are familiar with proper operation procedures and safety protocols outlined by the manufacturer.

From:

<https://bpm-wiki.cnsi.ucsb.edu/dokuwiki/> - **NSF BioPACIFIC MIP Wiki**



Permanent link:

https://bpm-wiki.cnsi.ucsb.edu/dokuwiki/doku.php?id=led_plate&rev=1728575343

Last update: **2024/10/10 15:49**